

User Manual

Version 20220110

Model: HGSM-XXX



Thank you for using the laser Cap. Please carefully refer to the user manual for safe and proper wearing instructions.


Intelligent Laser Hair Growth Cap

1. Brief Introduction to laser Hair Growth Cap

The laser cap is a soft laser hair growth product that not only helps control oil, prevent hair loss, and promote hair growth but also provides a gentle laser therapy experience that allows you to relax and enjoy 20 minutes of quiet time.

2. Technical Principle of the Intelligent Laser Hair Growth Cap

The Hair Growth Cap utilizes soft laser technology, also known as Low-Level Laser Therapy (LLLT). After applying LLLT, cellular metabolism increases, and the interaction between cells is enhanced. The soft laser, with its specific beam area, penetrates into the hair follicle layer, activating the hair follicle cells and improving their vitality.

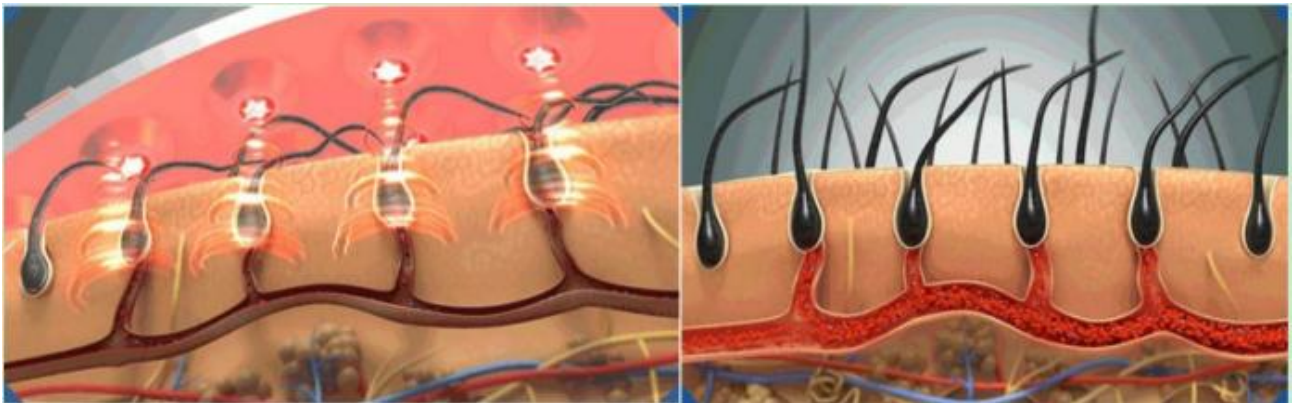

Stable and continuous low-energy laser irradiation stimulates capillaries in the head and increases blood flow

Increase the nutrient supply of hair follicles, renew hair follicle cells, and increase the toughness of hair

Before & After Use



Intelligent Laser Hair Growth Cap

3. Basic Product Information

Product list:

Item	P/N	product	Total beads pcs	Laser beads pcs	Red led pcs	Remark
1	HGSM-001	Laser Hair growth cap	306	306	0	
2	HGSM-002	Laser Hair growth cap	306	201	105	
3	HGSM-003	Laser Hair growth cap	82	82	0	
4	HGSM-004	Laser Hair growth cap	82	30	52	

Color: Black cap N.W: 250 grams G.W.: 532 grams

Packaging (L x W x H): 300 x 230 x 110 mm Output of laser diode: 680nm 5mw

Optional power bank: 5000mah




4. Usage Process:

- Check the power level of the power bank before use. If it has two or more bars of power, connect it to the hair growth cap.
- Put on the hair growth cap.
- The cycle time is fixed at 20 minutes/time; the cap will automatically power off when the time is up.
- Remove the hair growth cap.
- Neatly organize the hair growth cap and store it in the provided box. Charge the power bank for future using.

5. Contraindications:

- Do not use the hair growth cap if you are sensitive to light or sunlight due to your health condition or if you are taking photosensitive medication.

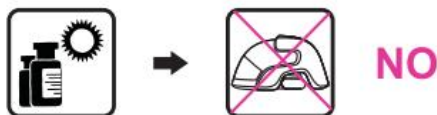

- Do not use this product if you have cancer on the scalp.



- Do not use this product if you have sensitive skin or are allergic to specific light at 680+/-10nm. Allergic skin or reactions may cause scalp itching or redness. If this occurs, please seek medical assistance.



NO

- Individuals under the age of 18 are not permitted to use this product.
- Individuals with hereditary hair loss or growth impairment are not recommended to use this product.

6 Warnings and precautions:

- Do not stare directly at the laser or use optical instruments such as magnifying glasses or microscopes to observe the laser.
- Stop using the product if you experience scalp itching or stinging during use.
- Do not wet the product or batteries.
- Keep the product away from high-temperature heat sources such as stoves, ovens, heaters, etc.
- Do not use when the hair or scalp is wet.
- Keep the product out of reach of minors and pets.
- Unplug the product before cleaning.
- Do not use the product while it is charging, please pay special attention.
- Do not attempt to repair the product yourself. In case of product malfunction, please contact customer service.
- Safely place the product to avoid dropping it, as dropping may cause malfunction or damage. If it falls and gets damaged, please contact customer service for assistance.

7. Cleaning and Storage

7.1 Cleaning and Maintenance

- Wipe the product with a soft cloth.
- Do not use corrosive chemical solvents to clean the product.
- Avoid liquid entering the product during cleaning to prevent functional failures.

7.2 Product Storage

After each use, please place the product in the dedicated packaging box. Store the packaging box in a safe environment, away from high temperatures, moisture, strong magnetic fields, and high acidity or alkalinity, to avoid damaging the product.

Intelligent Laser Hair Growth Cap

8. Package Contents

Item	Description	Qty	Unit	remark
1	Hair growth cap	1	pcs	
2	Power bank	1	pcs	Optional
3	Storage box	1	pcs	Non-disposable
4	User manual	1	pcs	



9. Warranty Card

Dear valued customer of products, thank you very much for choosing the hair growth cap and giving us the opportunity to provide you with the best service. To ensure our service meets your satisfaction, please carefully read this instruction and help keep this warranty card properly.

Warranty Terms:

1. During the warranty period, provides product coverage starting from the date of purchase and offers free repairs. Warranty repairs can only be performed at authorized service locations recognized by the brand. Please provide this warranty card or other valid proof of purchase with a purchase date when requesting warranty service.
2. Unauthorized disassembly or modification of the product is strictly prohibited. Violent operation of the product is also strictly prohibited. Unauthorized disassembly or violent operation will void the warranty.
3. Within the period starting from the date of purchase, the product will be repaired for any malfunctions that occur under normal usage conditions. However, the following circumstances, even within the warranty period, are not covered by the warranty and will be subject to paid repairs (component cost):
 - Product damage caused by violating the product's instructions for use.
 - Product damage caused by using uncertified extension components.
 - Product damage caused by exceeding the permissible usage environment, such as high temperatures or corrosive conditions.
 - Product damage caused by abnormal external forces, such as dropping or impact.
 - Counterfeit products.

Intelligent Laser Hair Growth Cap

- Inability to provide valid proof of purchase.
- Damages caused by force majeure events, such as fires.

To better serve you, please provide the following information:

Warranty Card	
Product Name	laser Cap
Product Model	
Customer Name	
Phone Number	
Date of Purchase	
Purchased from:	
Description of Malfunction	

Wish you all the best and extend our greetings!